




ENDURANCE® P2601
— BUILT TO LAST! —
TILT TOP OVERBED TRAY



It's ESSENTIAL for Good Health!

PARTS

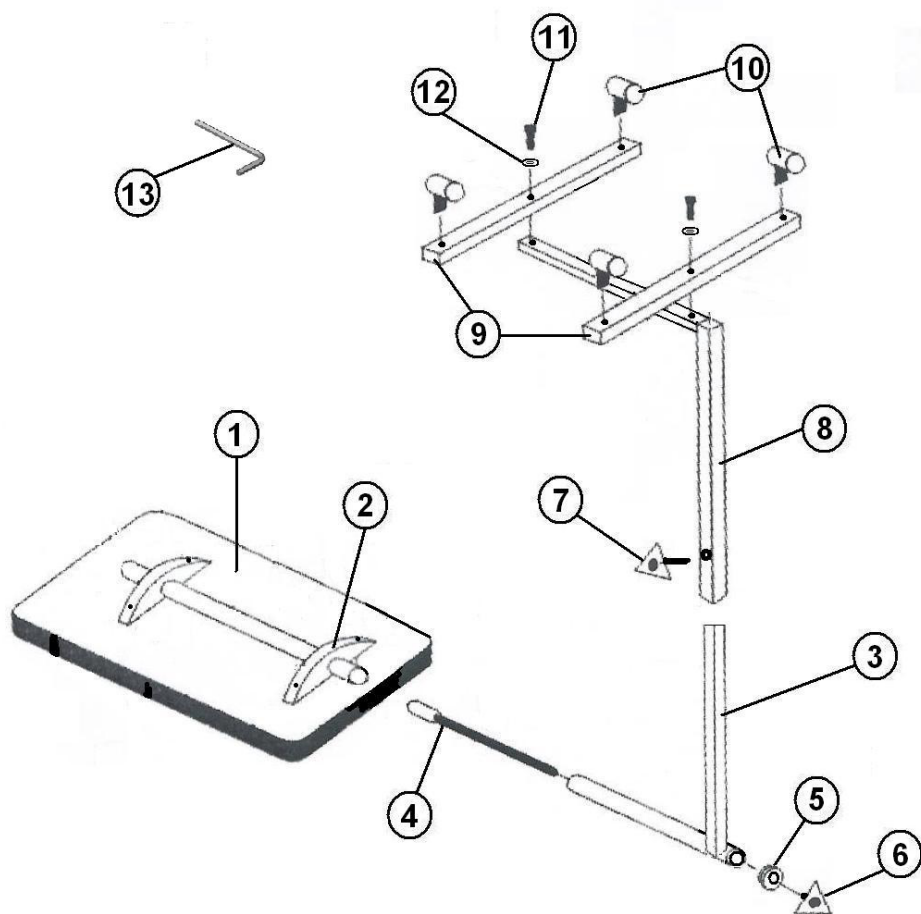

PARTS LIST

- | | |
|-----------------------------------|--------------------------------------|
| 1. Table Top x1 | 8. Table Support..... x1 |
| 2. Table Mount Bracket x1 | 9. Castor Support Arms x2 |
| 3. Table Mounting Arm x1 | 10. Castor with Lock x4 |
| 4. Angle Adjusting Screw x1 | 11. M8 x 1 1/4 x 34mm Socket |
| 5. Centering Bushing..... x1 | Head Cap Screw x2 |
| 6. Angle Adjusting Knob..... x1 | 12. 8 1/2 x 16 x 1mm Washer x2 |
| 7. Height Adjusting Knob x1 | 13. #6 Allen Wrench x1 |

P2601 OVERBED TABLE

ASSEMBLY INSTRUCTIONS

1. Invert Table Top (1). Slip the Angle Adjusting Screw Assembling (4, 5, 6) of the Table Mounting Arm (3) into the Table Mount Bracket (2) as shown and tighten the Angle Adjusting Knob (6) slightly.
2. Slip the arm of the Table Support (8) over the arm of the Table Mounting Arm (3) and secure it slightly with the Height Adjusting Knob (7).
3. Install Caster Support Arms (9) into horizontal arm of the Table Support (8), using M8 Socket Head Cap Screws (11) and Steel Washers (12) and #6 Allen Wrench (13).
4. Install all four Casters (10) into the Caster Support Arms (9) by pressing them into each hole of the Caster Support Arm (9) until secured in position.
5. Turn the Table assembly right side up and adjust as needed.
6. To adjust the Table Top (1) to desired height, lift the table and when in position, secure it with the Height Adjusting Knob (7).
7. To adjust the angle of the Table Top (10), loosen the Angle Adjusting Knob (6) and turn the Table Top(1) to desired position. Tighten the Angle Adjusting Knob (6) when done.

CAUTION

The Casters with Locks (10) shall be in locked position before using the Table.

Essential Medical Supply, Inc. Limited Warranty

Essential Medical Supply, Inc. ("Essential") warrants its products free from defects in materials and workmanship ("Warranty") for a period of three (3) years from the date of purchase ("Warranty Period"). The Warranty extends only to the original purchaser of Essential products. If within the Warranty Period there is a breach of Warranty as to an Essential product, such product shall be repaired or replaced at Essential's sole discretion. This Warranty does not include any labor or shipping costs incurred in replacement, part installation, or repair of any Essential product. **ESSENTIAL'S SOLE OBLIGATION AND YOUR EXCLUSIVE REMEDY UNDER THIS WARRANTY SHALL BE LIMITED TO REPAIR AND/OR REPLACEMENT OF THE DEFECTIVE ESSENTIAL PRODUCT.**

For Warranty service, please contact the dealer from whom you purchased your Essential product. In the event you do not receive satisfactory Warranty service, please write Essential directly at:

Essential Medical Supply, Inc.
Customer Service Department
6420 Hazeltine National Drive
Orlando, FL 32822-5121

Please provide the dealer's name, address, and date of purchase, and indicate the nature of the problem and the Essential product serial number.

Upon receipt of written notification, Essential will issue a return authorization. The defective unit or parts must be returned to Essential for Warranty inspection using the serial number as identification within thirty (30) days of the date of return authorization. **DO NOT** return products to Essential without Essential's prior written consent. Cash on delivery shipments will be refused and shipping charges must be prepaid.

Limitations and Exclusions: The Warranty shall not apply to products if 1) the Essential serial number has been removed or defaced; 2) the product has been damaged by negligence, accident, improper operation, maintenance or storage; and/or 3) the product has been modified without Essential's express written consent or repaired by an unauthorized third-party ("Warranty Disqualifying Events"). The Warranty Disqualifying Events include but are not limited to the addition of unauthorized parts or attachments or damage caused by circumstances beyond Essential's control, the occurrence of which shall be determined exclusively by Essential. The Warranty also shall not apply to problems or defects arising from normal wear and tear on Essential products and parts, including but not limited to rubber tips and hand grips, and failure to adhere to user instructions and safety precautions.

THE WARRANTY STATED HEREIN IS THE SOLE AND EXCLUSIVE WARRANTY OFFERED BY ESSENTIAL. THERE ARE NO OTHER WARRANTIES RESPECTING ESSENTIAL'S PRODUCTS, EITHER EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO ANY WARRANTY OF DESIGN, MERCHANTABILITY, OR FITNESS FOR A PARTICULAR PURPOSE, EVEN IF ESSENTIAL HAS BEEN INFORMED OF SUCH PURPOSE. NO AGENT OF ESSENTIAL IS AUTHORIZED TO ALTER OR EXCEED ESSENTIAL'S WARRANTY OBLIGATIONS AS SET FORTH HEREIN. ESSENTIAL SHALL NOT BE RESPONSIBLE FOR ANY INDIRECT, INCIDENTAL, SPECIAL, OR CONSEQUENTIAL DAMAGES ARISING OUT OF USE OF ESSENTIAL PRODUCTS.

Product Serial No. _____