

meywalk 2000

- Moves up and down during gait, allowing for a more normal heel strike, toe off and mid phase during gait
- Spring assist system to aide the user to “swing through” during gait
- **No more lifting to stand** the Meywalk has a “lever lift” system. The user transfers in at their wheelchair seat height and is brought into standing by pushing down on the lift bar



Wide range of accessories....



Midline leg separator



No Roll Back Brakes



Hip pads



Anti tippers (req'd with no back and drag brake options)



Drag Brakes

Other available options: castor swivel locks, height lowering kits

Finally... a Gait Trainer with Dynamic Movement to promote normal gait motion!

pacific rehab inc

Designs for function

Web www.pacificrehabinc.com
Phone 888-222-9040 480-595-6591
Fax 480-575-7907
Email information@pacificrehabinc.com
Address:
PO BOX 5406
Carefree, AZ 85377

Sizing	Small	Medium	Large
Ring Height	30-38"	33-40"	41-51"
Seat Height	19-26"	21-29"	27-37"
Max users weight	150 lb	220 lb	220 lb



Width of trunk req'd when ordering

Ease in Transfers!!!


To lower the seat and allow for transfer, first lock the unit, lower the tail gate and open the ring. Press the two silver buttons on the front of the walker and raise the lift bar which lowers the sit bar. After the user transfers, close the ring and tail gate, and push down on the lift bar.

The lever lift bar easily lifts the user into standing!



Dynamic Gait Movement

The Meywalk moves up and down with each step, allowing the user to come down into a heel strike and toe off. As the user steps forward, the spring assists with the swing through phase of gait

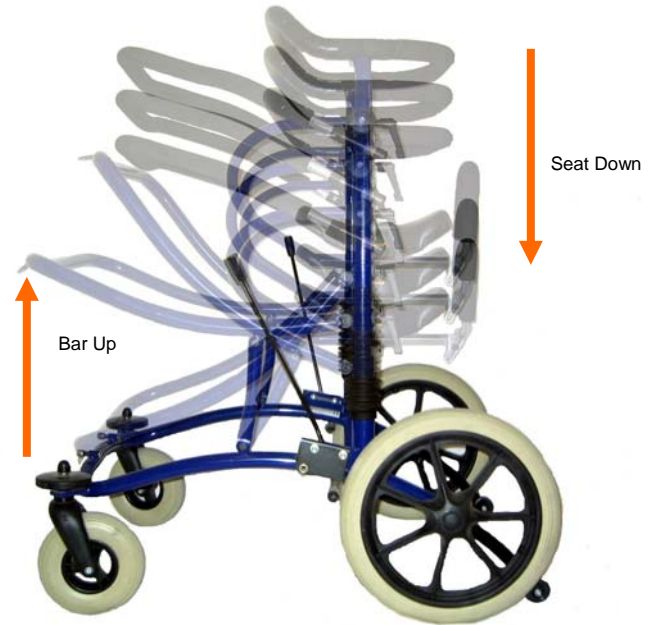

meywalk 2000

pacific rehab inc

Designs for function

Web www.pacificrehabinc.com
 Phone 888-222-9040 480-595-6591
 Fax 480-575-7907
 Email information@pacificrehabinc.com
 Address:
 PO Box 5406
 Carefree, AZ 85377